

CT/MR 2/3 Channel Endometrial Applicator Sets



Innovative technology and enhanced clinical advantages that facilitate CT and MR based brachytherapy planning.

Endometrium

Simplicity

The CT/MR 2/3 Channel Endometrial Applicator is designed for simplicity, even in the most advanced settings. It enables accurate and complete treatment with multiple lateral intrauterine tube spread distances and the ability to use a single applicator as either a 2 or 3 channel set. The single bracket connects all 3 intrauterine tubes with one lock knob and the key and lock design ensures proper lateral tube alignment.

Flexibility


The CT/MR 2/3 Channel Endometrial Applicator is available in 3 cm or 5 cm configurations, both of which can be used as either 2 or 3 channel applicators to effectively treat a variety of anatomies and disease states.

Reliability

The titanium construction with laser welded connections allows for enhanced stability and torque resistance. The key and lock design ensures proper lateral tube alignment in every case.

Features

- Autoclavable components
- Expected lifetime is 5 years
- Configurations for 3 cm and 5 cm spread
- 3.2 mm tube diameter
- Substantial dose coverage with a single set enabling two or three intrauterine tube configurations
- Simple, key lock bracket connection
- Applicator available with connections for BEBIG Medical afterloaders, Varian and Elekta/Nucletron transfer tubes*

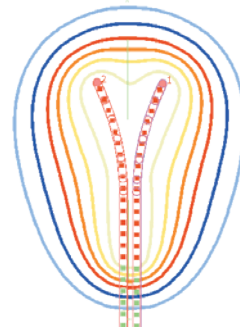

Isodose for 2-channel configuration



Sterilization box with select 2/3 Channel Endometrial Applicator Combo Set components

CT/MR 2/3 Channel Endometrial Applicator Sets

KBE-M1405

CT/MR 2/3 Channel Endometrial Applicator Combo Set for BEBIG Medical afterloaders

(Item numbers are included for reordering purposes)

CT/MR 2/3 Channel Endometrial Applicator Intrauterine Tube, right, L = 3 cm (KBE-1405-R-3)

CT/MR 2/3 Channel Endometrial Applicator Intrauterine Tube, left, L = 3 cm (KBE-1405-L-3)

CT/MR 2/3 Channel Endometrial Applicator Intrauterine Tube, central, L = 3 cm (KBE-1405-IU-3)

CT/MR 2/3 Channel Endometrial Applicator Intrauterine Tube, right, L = 5 cm (KBE-1405-R-5)

CT/MR 2/3 Channel Endometrial Applicator Intrauterine Tube, left, L = 5 cm (KBE-1405-L-5)

CT/MR 2/3 Channel Endometrial Applicator Intrauterine Tube, central, L = 5 cm (KBE-1405-IU-5)

Sealing Caps, (10 pcs each) (LAZ17-03)

Stainless Steel Ruler, L = 15 cm (KBE-7610)

Sterilization Box, 381 x 254 x 38 mm** (KBE-96-1522A)

Required Items

KBE-M1405, KBE-M1405-3, KBE-M1405-05

Qty.	Item	Part number
1	Transfer Tube, D = 3 mm, L = 1000 mm, Color code: green, Channel: 1	KBE-LAF1000-1
1	Transfer Tube, D = 3 mm, L = 1000 mm, Color code: yellow, Channel: 2	KBE-LAF1000-2
1	Transfer Tube, D = 3 mm, L = 1000 mm, Color code: blue, Channel: 3	KBE-LAF1000-3

KBE-M1405-5

CT/MR 2/3 Channel Endometrial Applicator 5 cm Spread Set for BEBIG Medical afterloaders

(Item numbers are included for reordering purposes)

CT/MR 2/3 Channel Endometrial Applicator Intrauterine Tube, right, L = 5 cm (KBE-1405-R-5)

CT/MR 2/3 Channel Endometrial Applicator Intrauterine Tube, left, L = 5 cm (KBE-1405-L-5)

CT/MR 2/3 Channel Endometrial Applicator Intrauterine Tube, central, L = 5 cm (KBE-1405-IU-5)

Sealing Caps, (10 pcs each) (KBE-17-03)

Stainless Steel Ruler, L = 15 cm (KBE-7610)

Sterilization Box, 381 x 254 x 38 mm** (KBE-96-1522A)

KBE-M1405-3

CT/MR 2/3 Channel Endometrial Applicator 3 cm Spread Set for BEBIG Medical afterloaders

(Item numbers are included for reordering purposes)

CT/MR 2/3 Channel Endometrial Applicator Intrauterine Tube, right, L = 3 cm (KBE-1405-R-3)

CT/MR 2/3 Channel Endometrial Applicator Intrauterine Tube, left, L = 3 cm (KBE-1405-L-3)

CT/MR 2/3 Channel Endometrial Applicator Intrauterine Tube, central, L = 3 cm (KBE-1405-IU-3)

Sealing Caps, (10 pcs each) (KBE-17-03)

Stainless Steel Ruler, L = 15 cm (KBE-7610)

Sterilization Box, 381 x 254 x 38 mm** (KBE-96-1522A)

Manufacturer:

Mick Radio-Nuclear Instruments, Inc.
 A BEBIG Medical Company
 521 Homestead Avenue
 Mount Vernon, NY 10550
 USA
 Phone 914 667 3999
 Fax 914 665 8834
 sales@micknuclear.com

Regional sales, marketing, and service in US

Klarity
 600 Industrial Parkway
 Heath, OH 43056
 USA
 Phone 740 788 8107
 Fax 740 788 8109
 info@klaritymedical.com
 www.klaritymedical.com



CT/MR 2/3 Channel Endometrial Applicator Sets

KBE-M0917

Set for Varian transfer tubes*

(Item numbers are included for reordering purposes)

CT/MR 2/3 Channel Endometrial Applicator Intrauterine Tube, right, L = 3 cm (KBE-0917-03-R3)

CT/MR 2/3 Channel Endometrial Applicator Intrauterine Tube, left, L = 3 cm (KBE-0917-03-L3)

CT/MR 2/3 Channel Endometrial Applicator Intrauterine Tube, central, L = 3 cm (KBE-0917-04-3)

CT/MR 2/3 Channel Endometrial Applicator Intrauterine Tube, right, L = 5 cm (KBE-0917-03-R5)

CT/MR 2/3 Channel Endometrial Applicator Intrauterine Tube, left, L = 5 cm (KBE-0917-03-L5)

CT/MR 2/3 Channel Endometrial Applicator Intrauterine Tube, central, L = 5 cm (KBE-0917-04-5)

Sterilization Caps, (10 pcs) (KBE-SRC093-12)

Stainless Steel Ruler, L = 15 cm (KBE-7610)

Sterilization Box, 381 x 254 x 38 mm** (KBE-96-1522A)

Required Items

3 Transfer Tubes. Please reach out to your Varian distributor for purchase.

KBE-M0916

Set for Elekta/Nucletron transfer tubes*

(Item numbers are included for reordering purposes)

CT/MR 2/3 Channel Endometrial Applicator Intrauterine Tube, right, L = 3 cm (KBE-0916-03-R3)

CT/MR 2/3 Channel Endometrial Applicator Intrauterine Tube, left, L = 3 cm (KBE-0916-03-L3)

CT/MR 2/3 Channel Endometrial Applicator Intrauterine Tube, central, L = 3 cm (KBE-0916-04-3)

CT/MR 2/3 Channel Endometrial Applicator Intrauterine Tube, right, L = 5 cm (KBE-0916-03-R5)

CT/MR 2/3 Channel Endometrial Applicator Intrauterine Tube, left, L = 5 cm (KBE-0916-03-L5)

CT/MR 2/3 Channel Endometrial Applicator Intrauterine Tube, central, L = 5 cm (KBE-0916-04-5)

Sterilization Caps, (10 pcs) (KBE-SRC093-12)

Stainless Steel Ruler, L = 15 cm (KBE-7610)

Sterilization Box, 381 x 254 x 38 mm** (KBE-96-1522A)

Required Items

3 Transfer Tubes. Please reach out to your Elekta/Nucletron distributor for purchase.

* The applicator has been tested for mechanical compatibility with the indicated manufacturer's transfer tubes to confirm there are no internal source obstructions.

Users are responsible for validating applicator compatibility, including treatment planning systems, with their afterloading equipment.

Varian is a trade name/trademark of Varian Medical Systems, Inc. Elekta/Nucletron are trade names/trademarks of Elekta AB (publ).

** Manufactured by PST Plastic Sterilizing Trays Corporation. Mick is a registered trademark of BEBIG Medical and its subsidiaries.

The mentioned products are not available in all markets. Please contact your local BEBIG Medical representative for more information.

Manufacturer:

Mick Radio-Nuclear Instruments, Inc.
A BEBIG Medical Company
521 Homestead Avenue
Mount Vernon, NY 10550
USA
Phone 914 667 3999
Fax 914 665 8834
sales@micknuclear.com

Regional sales, marketing, and service in US

Klarity
600 Industrial Parkway
Heath, OH 43056
USA
Phone 740 788 8107
Fax 740 788 8109
info@klaritymedical.com
www.klaritymedical.com

