


# CT/MR Contour Prostate Template Sets

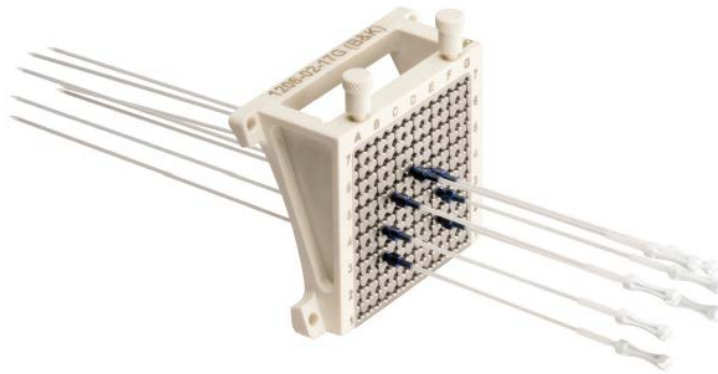

Facilitates interstitial needles for the treatment of prostate cancer.

## Prostate

### Simplicity

The CT/MR Contour Prostate Template facilitates the placement of interstitial needles for the treatment of prostate cancer. All templates are laser engraved with letters to denote the different vertical columns of needle channels and with numbers to denote the different horizontal rows of needle channels. The numbering/lettering layout corresponds to the type of ultrasound system being used.

### Flexibility

The CT/MR Contour Prostate Template is available to accommodate needles of different gauges, which are securely held in place by color-coded needle collets associated with the appropriate needle size. Both small and large template sizes are provided. The CT/MR Contour Prostate Template is curved to fit the sloping anatomy of the perineum. There are four suture holes provided for optional fixation of the template to the perineum.

### Reliability

The CT/MR Contour Prostate Template is constructed of PPSU (Polyphenylsulfone) and the template inserts and needle collets are constructed of titanium for continuous use and steam sterilization.

### Features

- Contoured to fit the sloping perineum
- Expected lifetime is 5 years
- MR conditionality allows for 1.5T or 3T
- Individual locking collet design enhances discrete needle targeting
- Available for variety of needle sizes: 15.5G, 17G, 18G, 1.65 mm, 2 mm
- Provided with four suture holes allowing the template to be directly sutured to the perineum
- Resistance-free bore holes for smooth needle entry
- Durable construction


Isodose




Template Needle Collets



Disposable Steel Needle Set

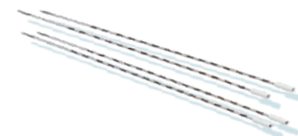

Plastic Needle Set

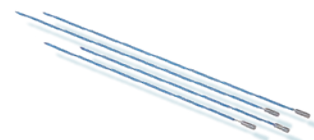

# CT/MR Contour Prostate Template Sets

## KBE-1204-02-XX (Small Matrix), for CIVCO Steppers, BK labeling scheme

(Item numbers are included for reordering purposes)

CT/MR Contour Prostate Template PPSU Body, small matrix, 17G  
CT/MR Contour Prostate Template Titanium Insert, small matrix, 17G  
Titanium Collet Wrench (KBE-1205-02)  
Mounting Pegs, (2 pcs) (KBE-1205-03)  
Mounting Bracket Lock Knobs, (2 pcs) (KBE-1205-12)  
Template Needle Collets (Titanium), seafoam, 17G,  
(6 pcs) (KBE-1105-01-EZ-17)

## KBE-1204-02-EZ-1.65 (Small Matrix), for CIVCO Steppers, BK labeling scheme

(Item numbers are included for reordering purposes)

CT/MR Contour Prostate Template PPSU Body, small matrix, 1.65 mm  
CT/MR Contour Prostate Template Titanium Insert, small matrix, 1.65 mm  
Titanium Collet Wrench (KBE-1205-02)  
Mounting Pegs, (2 pcs) (KBE-1205-03)  
Mounting Bracket Lock Knobs, (2 pcs) (KBE-1205-12)  
Template Needle Collets (Titanium), bronze, 1.65 mm,  
(6 pcs) (KBE-1105-01-EZ-1.65)

## KBE-1204-02-2.0 (Small Matrix), for CIVCO Steppers, BK labeling scheme

(Item numbers are included for reordering purposes)

CT/MR Contour Prostate Template PPSU Body, small matrix, 2.0 mm  
CT/MR Contour Prostate Template Titanium Insert, small matrix, 2.0 mm  
Titanium Collet Wrench (KBE-1205-02)  
Mounting Pegs, (2 pcs) (KBE-1205-03)  
Mounting Bracket Lock Knobs, (2 pcs) (KBE-1205-12)  
Template Needle Collets (Titanium), green, 2.0 mm,  
(6 pcs) (KBE-1105-01-2.0)

## KBE-1204-04-EZ-17 (Small Matrix), for CIVCO Steppers, Siemens labeling scheme

(Item numbers are included for reordering purposes)

CT/MR Contour Prostate Template PPSU Body, small matrix, 17G  
CT/MR Contour Prostate Template Titanium Insert, small matrix, 17G  
Titanium Collet Wrench (KBE-1205-02)  
Mounting Pegs, (2 pcs) (KBE-1205-03)  
Mounting Bracket Lock Knobs, (2 pcs) (KBE-1205-12)  
Template Needle Collets (Titanium), seafoam, 17G,  
(6 pcs) (KBE-1105-01-EZ-17)

## Optional Item

Item	Part number
Basket with lid for cleaning and sterilizing of small parts, 121 x 73 x 36 mm5	(KBE-07-01)

## KBE-1204-04-EZ-1.65 (Small Matrix), for CIVCO Steppers, Siemens labeling scheme

(Item numbers are included for reordering purposes)

CT/MR Contour Prostate Template PPSU Body, small matrix, 1.65 mm  
CT/MR Contour Prostate Template Titanium Insert, small matrix, 1.65 mm  
Titanium Collet Wrench (KBE-1205-02)  
Mounting Pegs, (2 pcs) (KBE-1205-03)  
Mounting Bracket Lock Knobs, (2 pcs) (KBE-1205-12)  
Template Needle Collets (Titanium), bronze, 1.65 mm, (6 pcs)  
(6 pcs) (KBE-1105-01- EZ-1.65)

## KBE-1204-04-2.0 (Small Matrix), for CIVCO Steppers, Siemens labeling scheme

(Item numbers are included for reordering purposes)

CT/MR Contour Prostate Template PPSU Body, small matrix, , 2.0 mm  
CT/MR Contour Prostate Template Titanium Insert, small matrix, 2.0 mm  
Titanium Collet Wrench (KBE-1205-02)  
Mounting Pegs, (2 pcs) (KBE-1205-03)  
Mounting Bracket Lock Knobs, (2 pcs) (KBE-1205-12)  
Template Needle Collets (Titanium), bronze, 2.0 mm,  
(6 pcs) (KBE-1105-01-2.0)

## KBE-1206-02-EZ-17 (Large Matrix), for CIVCO Steppers, BK labeling scheme

(Item numbers are included for reordering purposes)

CT/MR Contour Prostate Template PPSU Body, large matrix, 17G  
CT/MR Contour Prostate Template Titanium Insert, large matrix, 17G  
Titanium Collet Wrench (KBE-1205-02)  
Mounting Pegs, (2 pcs) (KBE-1205-03)  
Mounting Bracket Lock Knobs, (2 pcs) (KBE-1205-12)  
Template Needle Collets (Titanium), seafoam, 17G,  
(6 pcs) (KBE-1105-01-EZ-17)

### Manufacturer:

Mick Radio-Nuclear Instruments, Inc.  
A BEBIG Medical Company  
521 Homestead Avenue  
Mount Vernon, NY 10550  
USA  
Phone 914 667 3999  
Fax 914 665 8834  
sales@micknuclear.com

### Regional sales, marketing, and service in US

Klarity  
600 Industrial Parkway  
Heath, OH 43056  
USA  
Phone 740 788 8107  
Fax 740 788 8109  
info@klaritymedical.com  
www.klaritymedical.com



# CT/MR Contour Prostate Template Sets

## KBE-1206-02-EZ-1.65 (Large Matrix), for CIVCO Steppers, BK labeling scheme

(Item numbers are included for reordering purposes)

CT/MR Contour Prostate Template PPSU Body, large matrix, 1.65 mm  
 CT/MR Contour Prostate Template Titanium Insert, large matrix, 1.65 mm  
 Titanium Collet Wrench (KBE-1205-02)  
 Mounting Pegs, (2 pcs) (KBE-1205-03)  
 Mounting Bracket Lock Knobs, (2 pcs) (KBE-1205-12)  
 Template Needle Collets (Titanium), bronze, 1.65 mm, (6 pcs)  
 (KBE-1105-01-EZ-1.65)

## KBE-1206-02-2.0 (Large Matrix), for CIVCO Steppers, BK labeling scheme

(Item numbers are included for reordering purposes)

CT/MR Contour Prostate Template PPSU Body, large matrix, 2.0 mm  
 CT/MR Contour Prostate Template Titanium Insert, large matrix, 2.0 mm  
 Titanium Collet Wrench (KBE-1205-02)  
 Mounting Pegs, (2 pcs) (KBE-1205-03)  
 Mounting Bracket Lock Knobs, (2 pcs) (KBE-1205-12)  
 Template Needle Collets (Titanium), green, 2.0 mm, (6 pcs)  
 (KBE-1105-01-2.0)

## KBE-1206-04-EZ-17 (Large Matrix), for CIVCO Steppers, Siemens labeling scheme

(Item numbers are included for reordering purposes)

CT/MR Contour Prostate Template PPSU Body, large matrix, 17G  
 CT/MR Contour Prostate Template Titanium Insert, large matrix, 17G  
 Titanium Collet Wrench (KBE-1205-02)  
 Mounting Pegs, (2 pcs) (KBE-1205-03)  
 Mounting Bracket Lock Knobs, (2 pcs) (KBE-1205-12)  
 Template Needle Collets (Titanium), seafoam, 17G, (6 pcs)  
 (KBE-1105-01-EZ-17)

## KBE-1206-04-EZ-1.65 (Large Matrix), for CIVCO Steppers, Siemens labeling scheme

(Item numbers are included for reordering purposes)

CT/MR Contour Prostate Template PPSU Body, large matrix, 1.65 mm  
 CT/MR Contour Prostate Template Titanium Insert, large matrix, 1.65 mm  
 Titanium Collet Wrench (KBE-1205-02)  
 Mounting Pegs, (2 pcs) (KBE-1205-03)  
 Mounting Bracket Lock Knobs, (2 pcs) (KBE-1205-12)  
 Template Needle Collets (Titanium), bronze, 1.65 mm, (6 pcs)  
 (KBE-1105-01-EZ-1.65)

## KBE-1206-04-2.0 (Large Matrix), for CIVCO Steppers, Siemens labeling scheme

(Item numbers are included for reordering purposes)

CT/MR Contour Prostate Template PPSU Body, large matrix, 2.0 mm  
 CT/MR Contour Prostate Template Titanium Insert, large matrix, 2.0 mm  
 Titanium Collet Wrench (KBE-1205-02)  
 Mounting Pegs, (2 pcs) (KBE-1205-03)  
 Mounting Bracket Lock Knobs, (2 pcs) (KBE-1205-12)  
 Template Needle Collets (Titanium), green, 2.0 mm, (6 pcs)  
 (KBE-1105-01-2.0)  
 Sterilization Box, 381 x 254 x 38 mm1

### Required Items

KBE-1204-02-EZ-17, KBE-1204-04-EZ-17,  
 KBE-1206-02-EZ-17, KBE-1206-04-EZ-17

Qty.	Item	Part number
1	Disposable Steel Needle Set	KBE-SET0052
20	Transfer Tube, D = 3 mm, L = 1000 mm	KBE-LAF1000

### Required Items

KBE-1204-02-EZ-1.65, KBE-1204-04-EZ-1.65,  
 KBE-1206-02-EZ-1.65, KBE-1206-04-EZ-1.65

Qty.	Item	Part number
1	Plastic Needle Set	KBE-SET0012
20	Easy Click Transfer Tube, D = 3 mm, L = 1000 mm	KBE-LAG1000

### Required Items

KBE-1204-02-EZ-1.65, KBE-1204-04-EZ-1.65,  
 KBE-1206-02-EZ-1.65, KBE-1206-04-EZ-1.65

Qty.	Item	Part number
1	Plastic Needle, sterile, D = 2.0 mm, L = 250 mm, 15G3	KBE-LLA250-K15
5	CT/MR Obturator for 15G Needle, sterile, L = 250 mm3	KBE-LLH04-250B15
20	Easy Click Transfer Tube 6F, D = 3.0 mm, L = 1000 mm	KBE-LAG1000-6F

1 Manufactured by PST Plastic Sterilizing Trays Corporation  
 3 Manufactured by BEST Medical International, Inc.  
 5 Manufactured by Kögel GmbH

#### Manufacturer:

Mick Radio-Nuclear Instruments, Inc.  
 A BEBIG Medical Company  
 521 Homestead Avenue  
 Mount Vernon, NY 10550  
 USA  
 Phone 914 667 3999  
 Fax 914 665 8834  
 sales@micknuclear.com

#### Regional sales, marketing, and service in US

Klarity  
 600 Industrial Parkway  
 Heath, OH 43056  
 USA  
 Phone 740 788 8107  
 Fax 740 788 8109  
 info@klaritymedical.com  
 www.klaritymedical.com

